



Prof.ssa Elisa Ficarra
Prof.ssa Santa Di Cataldo
Eng. Marta Lovino
Eng. Alessio Mascolini

Politecnico di Torino
DAUIN
Dept. of Control and Computer Engineering

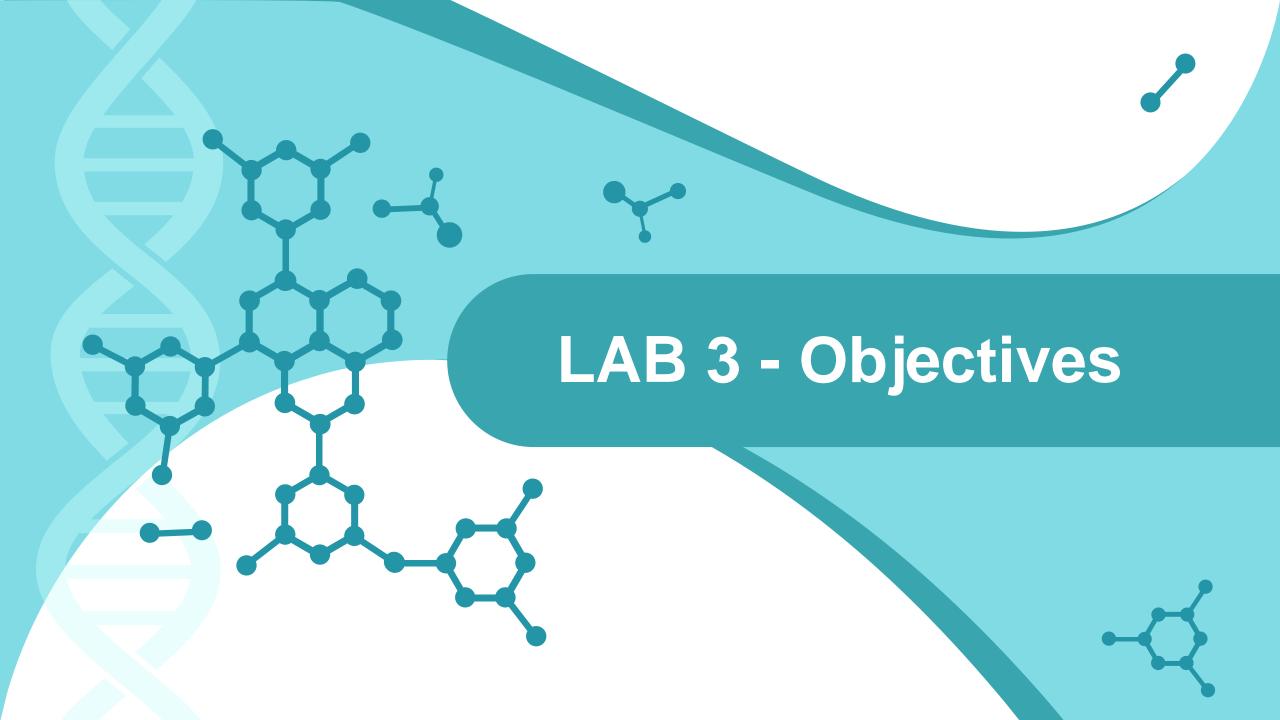

Objectives

- Understanding and manipulating different file formats (FASTA, FASTQ, SAM, BAM, GTF, VCF)
- Flags in SAM files
- Samtools, BCFtools examples

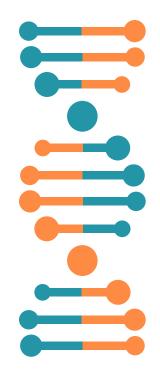

File formats – FASTA & FASTQ


FASTA

A sequence record in a FASTA format consists of a single-line description (sequence name), followed by line(s) of sequence data. The first character of the description line is a greater-than (">") symbol.

```
>read_id_0
GGTATGCTTCTGGGGCGGCAGTCGATAGGGCTAGACTCAGGTCCCGTGGC
>read_id_1
CACTGTGGCCCTCTTGGGGGGTGTCCACACGCCGCCCGTCGGCCCCCTCC
>read_id_2
GTTCTGTGGGTACCTCGCGGTTATGGTGTCGGGGGTATCCAAGGCACCCC
```



Similar to FASTA file, but with mapping quality information for each base. Both the sequence letter and quality score are each encoded with a single ASCII character.

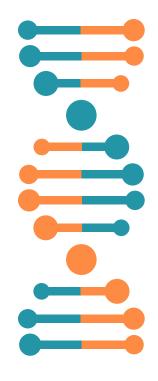

File formats – SAM

SAM stands for Sequence Alignment/Map format and it is the most common file output format for aligners. It is a TAB-delimited text format consisting of a header section, which is optional, and an alignment section.


```
@SQ
     SN:10
           LN:133797422
@SQ
     SN:18
           LN:80373285
     ID:bwa PN:bwa VN:0.7.17-r1188 CL:bwa mem
/home/marta/Documents/BIOINFORMATICS/BioInfoCourse/LAB alignment/tools/bwa index/bwa index
mate 1.fq mate 2.fq
ERX288614.1
                 10
                       55667614
                                        12S79M1I9M
                                                          55667850
     TTTTTCTAGACGGCAGGTCAGGTCCACCACTGACACGTTGGCAGTGGGGACACGGAAGGCCATGCCAGTGAGCTTCCCGTTCAGCTCAGG
MC:Z:19S82M
                                                     AS:i:41 XS:i:51
XA:Z:10,+15093507,101M,10;
ERX288614.1
           147
                 10
                       55667850
                                   39
                                        19582M =
                                                    55667614
318
     AGTCCTTCCACGATACCAAAGTTGTCATGGATGTCCTTGGCCAGGGGTGCTAAGCAGTTGGTGGTGCAGGAGGCATTGCTGATGATCTTG
B0<BB0<<<BB<B<<<< NM:i:4 MD:Z:14A39T8C0A17
                                    MC:Z:12S79M1I9M AS:i:62XS:i:53 XA:Z:10,-
91667278,15S83M3S,6;
```



File formats – SAM

- Remember that an aligner can report multiple alignments for the same read!! Depending from the application, this could be an issue to be fixed.
- Take a look here https://samtools.github.io/hts-specs/SAMv1.pdf for all details about this file format, paying particular attention to section 1.4 about mandatory fields (page 6).
- N.B. Genome positions in SAM files are in 1-based coordinate system

Col	Field	Type	Regexp/Range	Brief description
1	QNAME	String	[!-?A-~]{1,254}	Query template NAME
2	FLAG	Int	$[0, 2^{16} - 1]$	bitwise FLAG
3	RNAME	String	* [:rname:^*=][:rname:]*	Reference sequence NAME ⁹
4	POS	Int	$[0, 2^{31} - 1]$	1-based leftmost mapping POSition
5	MAPQ	Int	$[0, 2^8 - 1]$	MAPping Quality
6	CIGAR	String	* ([0-9]+[MIDNSHPX=])+	CIGAR string
7	RNEXT	String	* = [:rname:^*=][:rname:]*	Reference name of the mate/next read
8	PNEXT	Int	$[0, 2^{31} - 1]$	Position of the mate/next read
9	TLEN	Int	$[-2^{31}+1, 2^{31}-1]$	observed Template LENgth
10	SEQ	String	* [A-Za-z=.]+	segment SEQuence
11	QUAL	String	[!-~]+	ASCII of Phred-scaled base QUALity+33

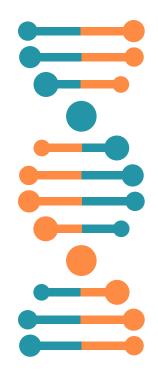

Flags in SAM file

Flags are used to keep track of alignment information in a compact way and uses 12 bits.

FLAG: Combination of bitwise FLAGs.⁷ Each bit is explained in the following table:

rom rigth		I	3it	Description		
	0	1	0x1	template having multiple segments in sequencing		
	1	2	0x2	each segment properly aligned according to the aligner		
	2	4	0x4	segment unmapped		
	3	8	0x8	next segment in the template unmapped		
	4	16	0x10	SEQ being reverse complemented		
	5	32	0x20	SEQ of the next segment in the template being reverse complemented		
	6	64	0x40	the first segment in the template		
	7	128	0x80	the last segment in the template		
	8	256	0x100	secondary alignment		
	9	512	0x200	not passing filters, such as platform/vendor quality controls		
	10	1024	0x400	PCR or optical duplicate		
	11	2048	0x800	supplementary alignment		



00000000001 --> 2°= 1 --> template having multiple segments in sequencing

 $0000000010 \longrightarrow 2^1=2 \longrightarrow each$ segment properly aligned according to the aligner

 $00000000100 --> 2^2 = 4 --> segment unmapped$

 $0001000000000 --> 2^{8} = 256 --> secondary alignment$

Bits can be combined:

 $00000001100 \longrightarrow 2^2 + 2^3 = 4 + 8 = 12 \longrightarrow segment unmapped and next segment in the template unmapped$



File formats – BAM

- BAM is the **binary version of a SAM file**. This means that BAM files are smaller than SAM files and this property is really helpful especially when we are working with huge files (e.g. a 30 GB SAM file can be compressed into a 17 GB BAM file). However, since BAM files are binary files, they are **not human readable**.
- N.B. Genome positions in BAM files are in 0-based coordinate system.



File formats – GTF

The Gene transfer format (GTF) is a file format used to **hold information about gene structure of a reference genome**. It is a tab-delimited text format based on the general feature format (GFF). This file format is really helpful when you want to know which biological feature (gene, exon,CDS, ...) is present in which genome positions. Genome positions in GTF files are in **1-based coordinate system**


```
gene id "ENSG00000223972";
        havana
gene_version "5"; gene_name "DDX11L1"; gene_source "havana"; gene_biotype
"transcribed_unprocessed_pseudogene"; havana_gene "OTTHUMG00000000961"; havana_gene_version "2";
        havana transcript
                                11869
                                        14409
                                                                         gene id "ENSG00000223972";
gene_version "5"; transcript_id "ENST00000456328"; transcript_version "2"; gene_name "DDX11L1";
gene source "havana"; gene biotype "transcribed unprocessed pseudogene"; havana gene
"OTTHUMG00000000961"; havana_gene_version "2"; transcript_name "DDX11L1-002"; transcript_source
"havana"; transcript biotype "processed transcript"; havana transcript "OTTHUMT00000362751";
havana_transcript_version "1"; tag "basic"; transcript_support_level "1";
                                                                 gene id "ENSG00000223972";
        havana <mark>exon</mark>
gene version "5"; transcript id "ENST00000456328"; transcript version "2"; exon number "1"; gene name
"DDX11L1"; gene_source "havana"; gene_biotype "transcribed_unprocessed_pseudogene"; havana_gene
"OTTHUMG00000000961"; havana_gene_version "2"; transcript_name "DDX11L1-002"; transcript_source
"havana"; transcript_biotype "processed_transcript"; havana_transcript "OTTHUMT00000362751";
havana_transcript_version "1"; exon_id "ENSE00002234944"; exon_version "1"; tag "basic";
transcript support level "1";
"1";
```



File formats – VCF

VCF is a text file format (most likely stored in a compressed manner) and it holds information about reads. It contains meta-information lines, a header line, and then data lines each containing information about a position in the genome. The format also has the ability to contain genotype information on samples for each position.

#CHROM	POS	/II		REF A	LT	QUAL	FILTER INFO FORMAT sorted.bam		
10 48	8636) <u>†</u>	<*>	0		DP=1;I16=0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0;QS=0,0;MQ0	F=0 PL	0,0,0
10 48	8637	•	G	<*>	0		DP=2;I16=0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0;QS=0,0;MQ0	F=0 PL	0,0,0
10 48	8638	•	Α	<*>	0	•	DP=2;I16=1,0,0,0,70,4900,0,0,60,3600,0,0,2,4,0,0;QS	=1,0;MQ0F=0	PL
0,3,60									
10 48	8639		Α	<*>	0		DP=2;I16=1,0,0,0,74,5476,0,0,60,3600,0,0,3,9,0,0;QS	=1,0;MQ0F=0	PL
0,3,60									
10 48	8640		G	<*>	0		DP=2;I16=1,0,0,0,70,4900,0,0,60,3600,0,0,4,16,0,0;Q	S=1,0;MQ0F=0	PL
0,3,60									
10 48	8641		Α	<*>	0		DP=2;I16=1,0,0,0,70,4900,0,0,60,3600,0,0,5,25,0,0;Q	S=1,0;MQ0F=0	PL
0,3,60									
10 48	8642	•	C	T,<*>	0		DP=2;I16=0,0,1,0,0,0,70,4900,0,0,60,3600,0,0,6,36	;QS=0,1,0;SG	B=-
3.37988	5;MQ0F=6	9	PL	60,3,0	,60,3,	60			
	10 4 10 4 0,3,60 10 4 0,3,60 10 4 0,3,60 10 4	10 48636 10 48637 10 48638 0,3,60 10 48639 0,3,60 10 48640 0,3,60 10 48641 0,3,60 10 48642	10 48636 . 10 48637 . 10 48638 . 0,3,60 10 48639 . 0,3,60 10 48640 . 0,3,60 10 48641 .	10 48636 . T 10 48637 . G 10 48638 . A 0,3,60 10 48639 . A 0,3,60 10 48640 . G 0,3,60 10 48641 . A 0,3,60 10 48642 . C	10 48636 . T <*> 10 48637 . G <*> 10 48638 . A <*> 10 48638 . A <*> 0,3,60 10 48639 . A <*> 0,3,60 10 48640 . G <*> 0,3,60 10 48641 . A <*> 0,3,60 10 48641 . A <*> 0,3,60 10 48642 . C T,<*>	10 48636 . T <*> 0 10 48637 . G <*> 0 10 48638 . A <*> 0 10 48638 . A <*> 0 0,3,60 10 48640 . G <*> 0 0,3,60 10 48641 . A <*> 0 0,3,60 10 48641 . A <*> 0 0,3,60 10 48642 . C T,<*> 0	10 48636 . T <*> 0 . 10 48637 . G <*> 0 . 10 48638 . A <*> 0 . 10 48638 . A <*> 0 . 10 48639 . A <*> 0 . 10 48640 . G <*> 0 . 10 48641 . A <*> 0 . 10 48641 . A <*> 0 . 10 48642 . C T,<*> 0 . 10 48642 . C T,<*	DP=1;I16=0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,	10 48636 . T <*> 0 . DP=1;I16=0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,0,

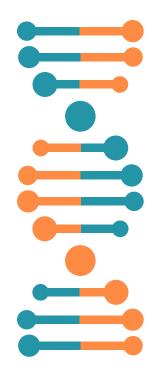


Tools and data

INSTALL SAMTOOLS, BCFTOOLS and download human gtf file

conda activate Bioinfo_labs (or source activate Bioinfo_labs)
conda install -c bioconda samtools
conda install -c bioconda bcftools

Download **Homo_sapiens.GRCh38.95.gtf.gz** file from ftf.gz, move to the correct folder and extract it:

gunzip -d Homo_sapiens.GRCh38.95.gtf.gz

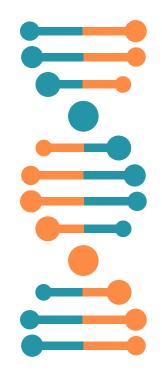

Samtools, really basic usage

SAM Tools provide various utilities for manipulating alignments in the SAM format, including sorting, merging, indexing and generating alignments in a per-position format.

SAM/BAM conversions


As previously pointed out, BAM format is the binary version of a SAM file. The conversion can be performed using samtools:

samtools view -S -b my.sam > my.bam

BAM sorting

When multiple selections have to be performed onto a huge file, it is convenient to sort that file according to certain criteria (e.g. genomic region) in order to search for the required information in a faster way. For BAM sorting you can use:

samtools sort my.bam > my-sorted.bam

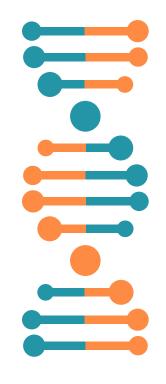

Samtools to filter SAM files using FLAGS

Among the many potentials of samtools view there is that of filtering the reads using the FLAG field in SAM files. To have an overall idea of how samtools view works open your terminal, activate Bioinfo_labs environment and type

samtools view

to get its manual page.


Take a look at -f and -F options. These options allow us to filter reads in a SAM file following a combination of criteria relying on alignment flags.

- -f INT only include reads with all of the FLAGs in INT present
- -F INT only include reads with none of the FLAGS in INT present

E.g. to obtain a SAM file with no unmapped reads and no secondary alignments we can exploit bit in position 2 (2^2 =4) and bit in position 8 (2^8 =256).

samtools view -F 260 bwa_out.sam > unique_aligned.sam

-F 260 $(2^2 + 2^8)$ means that in the final SAM file will be printed only reads for which bit number 2 or bit number 8 is not set to 1.

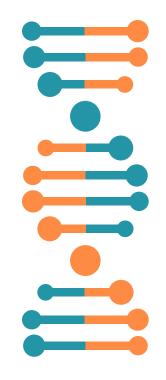


Bcftools


BCFtools is a set of utilities that manipulate variant calls in the Variant Call Format (VCF) and its binary counterpart BCF.

BCFTOOLS to create VCF file

Use bcftool to convert a sorted BAM file into VCF:

bcftools mpileup --fasta-ref reference_chr10_chr18.fa sorted.bam > sorted.vcf

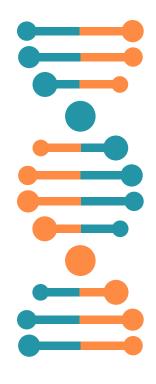


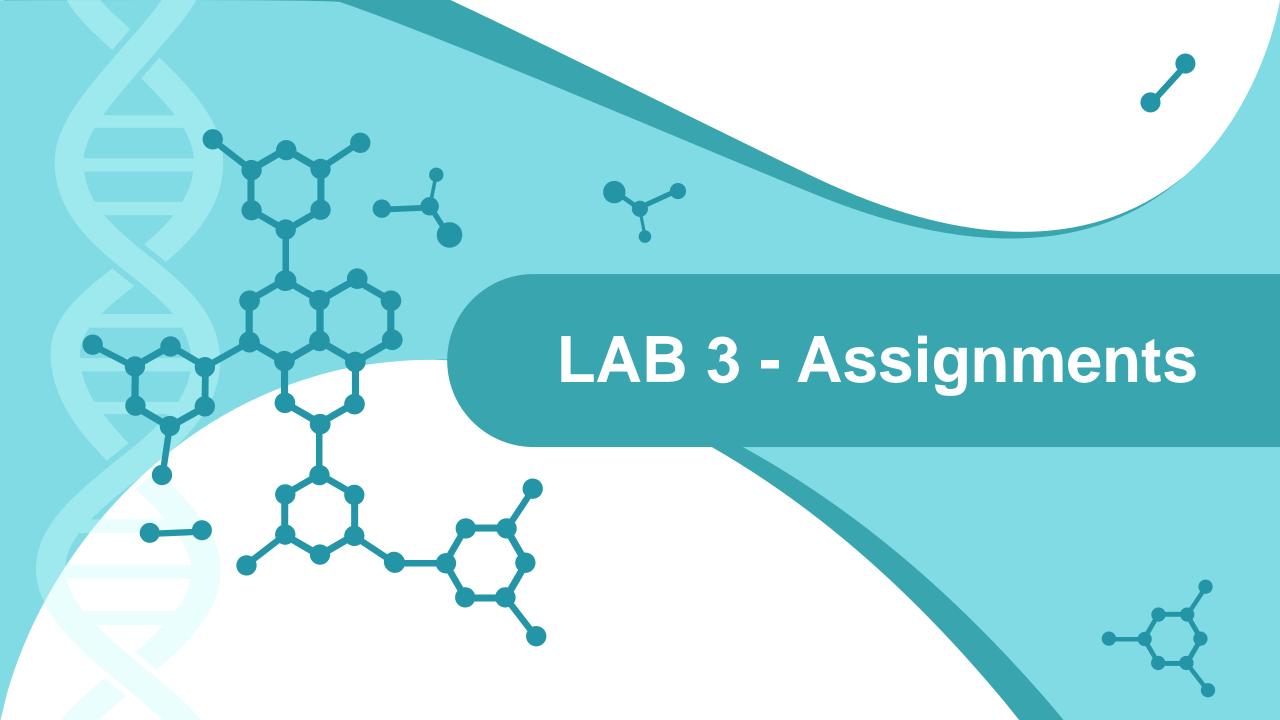

Assignment 1: Search for SNP and deletions in a sample

In order to obtain biological information on SNPs, insertions and deletions, manipulate SAM files in the following way:

- 1. Convert SAM file obtained in LAB2 to BAM
- 2. Sort BAM to make the following process faster.
- 3. Use beftools to obtain a VCF file
- 4. Write a Python program to parse the VCF file and obtain only Single Nucleotides Polymorphism (SNP) for which the information is complete.

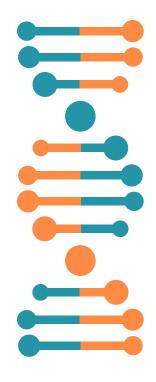


Assignment 1: Search for SNP and deletions in a sample

5. Write a Python program to parse VCF file and obtain only Insertions or Deletions (INDEL) for which the information is complete

```
e.g.

10 3643016 . GCTC G 0 .


INDEL; IDV=2; IMF=1; DP=2; I16=0,0,1,1,0,0,92,5800,0,0,14,98,0,0,50,1250; Q
S=0,1; VDB=0.06; SGB=-0.453602; MQSB=1; MQ0F=0 PL 14,6,0
```

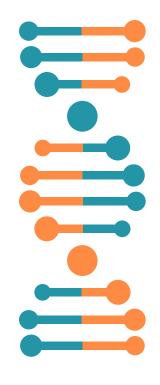

Assignment 2: Raw read count for protein coding genes

With this assignment we want to calculate how many reads of our SAM file have been mapped on the protein coding genes of chromosome 10 and 18. First of all, we have to obtain names and genomic positions of each protein coding gene from chr10 and chr18.

- Write a python program that parses the gtf file and extract information about chromosome 10 and 18 (means first column equal to 10 or 18). Select only rows for which the feature is equal to "gene" and gene_biotype equal to "protein_coding". Alternatively, you can perform this step using awk bash command (see awk documentation online https://www.gnu.org/software/gawk/manual/gawk.html#Very-Simple).
- Use samtools view to filter unmapped reads and supplementary alignments from the SAM file you obtained in LAB2. Save your results into a new SAM file (e.g. unique_aligned.sam)
- Using information from the reduced SAM and gtf files, write a Python script to calculate the raw read count for each protein coding gene in chr10 and chr18 (since the original files we provided you were truncated for computational feasibility, you will find a lot of genes with read count equal to zero).

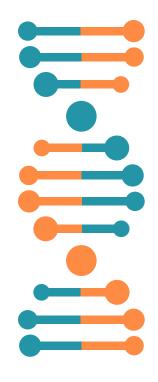

LAB3 – Take home message

- SAM files are usually from 5 to 30 GB each!! Read these files line by line
- BAM is the compressed version of a SAM file. The compression depends on the specific SAM file, but it can vary from 1/2 to 1/10.
- GTF is about 1 GB. You can store it entirely in the memory
- VCF file depends on the SAM file and on the number of alterations in the sample. Usually from 5 to 30 GB!!



Questions?

Remember: no question is stupid